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<u>CZ</u>	<u>3</u>
<u>SK</u>	<u>6</u>
<u>EN</u>	<u>9</u>

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NUV041310

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NUV041310



DEAR CUSTOMER

THANK YOU FOR CHOOSING THIS BABY SCALE. PLEASE READ THESE INSTRUCTIONS FOR USE CAREFULLY AND KEEP THEM FOR FUTURE USE.

SPECIFICATIONS

- Accurately measures weight in 5g increments up to 20kg
- HOLD feature captures precise weight measurement of an active baby
- TARE feature subtracts added weight of diaper, clothing or blanket
- Easy-to-read, white backlit LCD display
- Unique touch button design
- Includes 2xAAA batteries
- Auto power off and low battery indication

DISPLAY

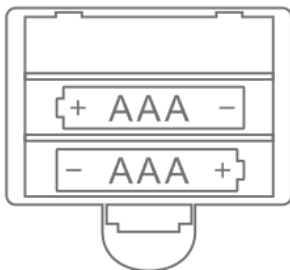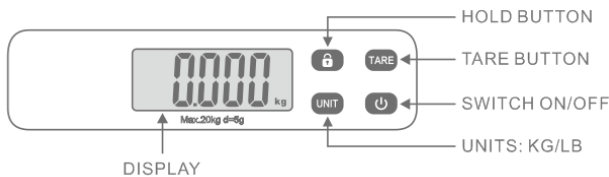POWER SUPPLY

Please insert batteries or remove the insulating foil from the batteries before using. When "Lo" appears in the display, the batteries need to be replaced.

OPERATING INSTRUCTIONS

The weight is displayed in steps of 5g from a minimum weight of 30g. Press the button "⏻" to switch on the scale. When "0.000kg" appears on the display, the scale is ready for measurement. Place the baby on the scale carefully. Take particular care with the baby's head.

HOLD FUNCTION:

The scale has a “**HOLD**” function: **AUTO HOLD, MECHANICAL HOLD.**

AUTO HOLD: When the baby is placed on platform, he/she may move around. When he/she is still enough, the AUTO HOLD function will establish the weight on display. The sign “” flashes and then weight is established.

MECHANICAL HOLD: If the baby is moving around too vigorously, you may select the MECHANICAL HOLD function by pressing “” button. The sign “” flashes and weight is established on display. The sign “” remains on display until you press “” button or “TARE” button. Then it shows “0.000kg” which means the scale is ready for next weighing.

TARE FUNCTION:


If you prefer to place a towel or blanket on the scale and weigh the baby, please do so BEFORE switching the scale on. If the towel or blanket weighs over 30g, please press TARE button.

“TARE” is shown on display. Then place the baby onto towel or blanket. It shows baby’s weight.

UNITS:

The scale is with units of kilograms(kg) and pounds(lb), which can be switched by pressing “UNIT” button.

SWITCH OFF

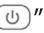After weighing is finished, you can switch off the scale by pressing “” button. Or, the scale will automatically switch off if no weighing is carried out within 180 seconds.

OVERLOAD INDICATION

The scale can accept a maximum load of 20 kg. When overloaded, “EEEE” sign will appear on display.

CARE AND ATTENTION

- Place the scale on a hard, flat surface to ensure greatest accuracy.
- Cleaning: You can clean the scale with a damp cloth and a little washing up liquid if necessary. Never immerse the scale in water or rinse it under running water.
- Protect the scale against knocks, damp, dust, chemicals, drastic changes in temperature, electromagnetic field and source of heat that are too near.
- Keep packaging material away from children.
- DO NOT dispose of the old batteries in fire. The batteries may explode or leak. Recycle or manage them as hazardous waste.
- If the scale is not used in a long period, please take the batteries out.

WARRANTY CERTIFICATE

1. The seller provides a product warranty for 24 months from the date of sale to the consumer.
2. The seller is obligated to complete the warranty certificate upon sale with the date and stamp of the store. It is not permitted to date and stamp the warranty certificate after the product is sold.
3. In the event of a complaint, the consumer must present the original Warranty Certificate, the original product purchase receipt and make a record of the claimed event (a description of the circumstances in which the defect occurred or a brief description of the defect with the consumer's request).
4. The product on the claim will only be accepted if it's clean.
5. The warranty period shall be extended when the consumer has filed a claim for warranty repairs until the time the consumer is obliged to take on the product after repairs are completed.
6. The warranty expires if the consumer or other person, after purchase, does not follow instructions during assembly, maintenance, cleaning and use in accordance with the instructions for use contained in these warranty terms. The consumer is not entitled to make any changes to the product.
7. For mechanical damage due to improper use, in particular scratches, breakage of plastic parts, cracked or broken joints, discolouration of cover colours resulting from non-compliant cleaning or brightening caused by exposure of the product to direct sunlight or moisture, a claim will not apply. It is also not possible to claim for cracked parts due to a fall, abrasion or tearing of the cover and wear-and-tear of the product by normal use. The claim cannot be applied if product parts are stolen, lost or damaged by their own fault. These defects will be repaired by the consumer, including shipping costs.
8. The buyer was shown how to assemble the car seat upon purchase or it was assembled by the seller. (Cross out what is not applicable) and the buyer was familiar with the instructions for use when buying the item. The provision of the guarantee is without prejudice to the buyer's rights, which are subject to the purchase of the item in accordance with special legal regulations.

Distributor: BabyDirekt s.r.o.
Masarykova 118, 664 42 Modřice, Czechia
e-mail: info@babydirekt.cz
www.babydirekt.eu